

Operator's Manual

Model **614B** • Laboratory Centrifuge

Table of Contents

| | |
|---|---------|
| Model Description | pg. 2 |
| Supplied Equipment | pg. 2 |
| Features | pg. 3 |
| Specifications | pg. 3 |
| Setup Location and Procedure | pg. 4 |
| Operation | pg. 5 |
| Rotor Removal and Installation | pg. 5 |
| Tube Holder Configurations | pg. 6 |
| Care and Preventative Maintenance | pg. 7,8 |
| Safety..... | pg. 8 |
| Calibration and Ground Testing | pg. 8 |
| Troubleshooting | pg. 8 |
| Replacement Parts | pg. 9 |
| Available Accessories | pg. 9 |
| Warranty..... | pg. 10 |

WARNING: For the safety of both the operator and service personnel, care should be taken when using this centrifuge if handling substances that are known to be toxic, radioactive or contaminated with pathogenic microorganisms. When Risk Group II materials are used, (as identified in the World Health Organization "Laboratory Bio-Safety Manual"), a Bio-Seal should be employed. The Bio-Seal accessory for the model 614 tube holders is the non-aerosol shield cap, p/n 7713011. In the event that materials of a higher risk group are being used, more than one level of protection must be provided. The use of flammable or explosive materials as well as those materials which have a vigorous chemical reaction is prohibited. For your safety and the durability of your machine, never transport or store centrifuge with tube holders inside machine.

NOTES:

Model Description:

The **Model 614B** is a continuous duty centrifuge designed for the small lab or doctor's office for the purposes of separating laboratory fluids. The unit is controlled by a mechanical timer settable from 1 to 30 minutes. Samples can be safely viewed through the transparent lid. In the event that the lid is opened during a run, the power to the motor is disconnected.

For warranty information, turn to page 10.

Supplied Equipment*:

The following items come standard with each **Model 614B** centrifuge:

- | | | |
|--|-----|---------|
| 1. One (1) six-place fixed-angle rotor | p/n | 7786047 |
| 2. Operator's manual | | |
| 3. Six (6) 125mm tube holders | p/n | 7713032 |
| 4. Six (6) 100mm tube holders | p/n | 7713031 |
| 5. Six (6) 75mm tube holders | p/n | 7713033 |



* The rotor and rotor accessories are rated for a rotation frequency of 3,500 RPM.

WARRANTY:
 The Drucker Company
 warranties that this centrifuge
 is free from defects in
 workmanship and parts
 for 12 months.

Features:

- Fixed-angle rotor for spinning 75mm, 100mm and 125mm test tubes in specially designed tube holders
- Cool-Flow air flow design that prevents overheating of samples
- Heavy gauge steel construction for safety and durability
- Lid safety switch that prevents the centrifuge from operating unless the lid is closed and latched
- Removable rotor for easy cleaning
- Brushless A/C motor
- Clear lid for safe observation of samples and optical calibration of speed

Specifications:

General Specifications for the Model 614B Centrifuge

| | |
|--|---|
| Nominal Speed (125mm holders): | 3,150 (± 100) RPM |
| Nominal RCF (125mm holders): | 1,200 (± 80) xg |
| Nominal Speed (100mm holders): | 3,250 (± 100) RPM |
| Nominal RCF (100mm holders): | 1,080 (± 80) xg |
| Nominal Speed (75mm holders): | 3,300 (± 100) RPM |
| Nominal RCF (75mm holders): | 950 (± 80) xg |
| Maximum capacity: | 90 mL (6 x 15 mL)* |
| Overall Dimensions (H x W x D): | 8.5 in. x 11 in. x 12.5 in. |
| Centrifuge Motor: | 1/30 HP, AC |
| Nominal Acceleration Time: | 45 seconds |
| Protection Breaker: | 4 Amp. re-settable |
| Timer: | mechanical, 1 to 30 minutes accuracy ± 10% |
| Current Requirement: | 1.0 Amps |
| Voltage Requirement: | 115 (±10) Volts |
| Frequency: | 60 Hz |
| Weight: | 18.5 lbs. |

Before using any cleaning or decontamination methods except those recommended by the manufacturer, users should check with the manufacturer that the proposed method will not damage the equipment. See page 7, (bottom), for the recommended cleaning solutions.

CE } ^Á ~ • ^Á [c@^/Ác@æ}Ác@ [• ^Á •] ^Áã, ^ÁãÁà ^Ác@^Á Tæ } ~æ&c ~!^/Áã • Á ^Áç] [ã&ac| ^Á] ; [@ã&ac^ÁãÉÁÁ

* Maximum sample density is 1.15 grams / mL,
 (water density = 1.0 grams / mL)



Tæã^Áã }Ác@^ÁWÜÇÉÁã

200 Shadylane • Philipsburg, PA 16866
 Phone: 814-692-7661 • Fax: 814-692-7662 • www.druckercompany.com

Setup Location:

1. Unpack the centrifuge and verify that all of the supplied equipment is present.
2. Choose a setup location which meets the following criteria:
 - a) A bench top clearance height of 20" is required in order to open the lid.
 - b) The clearance envelope is the space around the centrifuge which is required for safety. Choose a setup location which will allow for a clearance envelope of at least 24" x 24", (with the centrifuge at the center). No person or hazardous material shall be permitted in the clearance envelope during operation. The operator time within the envelope shall be limited to the time necessary for loading, unloading and centrifuge operation only.
 - c) Proper ventilation is necessary to prevent the overheating of samples as well as premature failure of the centrifuge. Choose an area which will allow unencumbered air flow.
 - d) The centrifuge is designed to secure to the operating surface by four suction feet. No adjustment is necessary for leveling the centrifuge, however, the surface should be flat and level.
 - e) **Be sure the outlet is always within reach as the line cord is the means of emergency disconnection!**

Initial Setup Procedure:

If any problems are found during the initial setup procedure, refer to the troubleshooting section on page 8. For further assistance, contact The Drucker Company at 814-692-7661.

1. Plug the centrifuge in to an approved electrical outlet. For electrical safety, the unit must always be properly grounded.
2. Turn the latch counter-clockwise and open the lid.
3. Spin the rotor by hand; check for free and level rotation.
4. Close the lid. Rotate the lid knob clockwise to its complete stop position.
5. Turn the centrifuge on by turning the timer to 10 minutes.
6. Listen to the centrifuge. A smooth whirring sound should be heard.

After the centrifuge has passed this procedure it is ready for operation.

BALANCED LOADS

Your centrifuge must contain a balanced load in order to work properly. Use the following rules when loading the rotor. Spinning balanced loads will extend the life of the machine and produce better results.

1. *Opposing tube holders must be identical and must contain the same cushion, or none at all.*
2. *Opposing tube holders must be empty or loaded with equally weighted samples.*


match the weight of the unpaired sample and place it across from this sample.

Replacement Parts:

| Part No. | Description |
|----------|-------------------------------------|
| 7724037 | Foot, rubber |
| 7751068 | Switch, lid safety |
| 7786047 | Rotor, six-place fixed-angle |
| 7735050 | Motor, 1/30 H.P. |
| 7722027 | Timer, mechanical |
| 7751043 | Circuit Breaker |
| 7724022 | Front Panel Label |
| 7760002 | Power cord |
| 7714101 | Pawl, latch, lid |
| 7714103 | Knob, latch, lid |
| 7712260 | Lid |
| 7724071 | Hinge, friction |
| 7732018 | Seal, lid gasket |
| 7732019 | Seal, rotor chamber gasket |
| 7713033 | Tube holder, green, for 75mm tubes |
| 7713031 | Tube holder, red, for 100mm tubes |
| 7713032 | Tube holder, black, for 125mm tubes |

Available Accessories:

1" Tube cushion

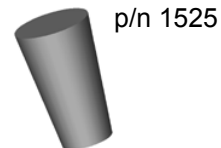

Shield caps
p/n
7713011



200 Shadylane • Philipsburg, PA 16866


Phone: 814-692-7661 • Fax: 814-692-7662 • www.druckercompany.com

Tube Holder Configurations:

The fixed-angle rotor that came with your **Model 614B** is capable of spinning test tubes up to 17mm x 125mm. Use the following chart and drawing to determine which tube holder and cushion combination should be used with your application.

DIRECTIONS:

1. Compare the tube to be spun with the four boxes shown below.
2. Find the box that most closely matches the tube's length. **NOTE:** The tube length with its stopper or cap must be shorter than the chosen box or the tube will not fit properly in the tube holder.
3. Match the letter from the chosen box with one of the configurations shown.

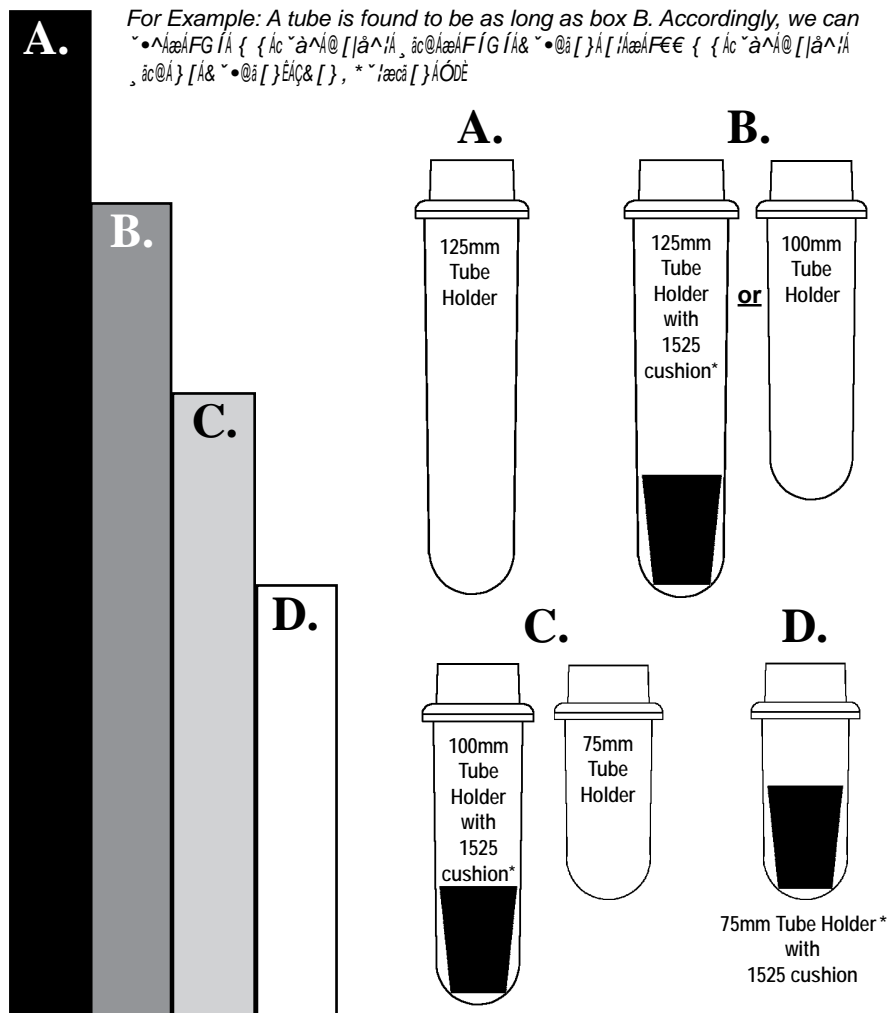

ÉÁ V@á•Á]æ /cá•kæçæi]æà/Áæ•kæ }kæ&&•• [/ ÉÁÁÖ []cæ&cá V@ÁÖ; ~\^ÁÖ [{]æ } ^- [/kæ••j•cæ }&^ÉÁ

Care and Preventative Maintenance:

With proper care and maintenance your centrifuge will provide years of laboratory service.

For proper care, the following steps should be taken:

1. **Provide Adequate Ventilation:** For cooling purposes, the centrifuge draws in ambient air through the air intake cover on the top of the lid and exhausts this air in the rear of the base. The centrifuge should be placed on a hard smooth surface for good air circulation.
2. **Always Spin Balanced Loads:** Make certain that you are always spinning a balanced load. The Model 614B has a unique counter balanced motor mounting design which, along with its rubber suction feet, produces excellent vibration dampening. However, out-of-balance loads may break glass test tubes and may produce unsatisfactory separation results. Proper load balancing will improve sample separation and extend the life of the centrifuge. Refer to page 4 on balanced loads for additional information on balancing the load.
3. **Keep the Tube Holders Clean:** **NOTE:** Always follow the safety guidelines of your laboratory to properly clean up and/or dispose of materials in the event that a substance known to be potentially toxic, radioactive or contaminated with a pathogenic microorganism is spilt in or on the centrifuge. Small glass fragments left in the tube holder after a tube breakage may adhere to the next test tube inserted in that holder. When this tube is handled, these fragments may puncture protective gloves and lacerate the operator's fingers or hand. Remaining fragments may provide stress points on subsequent tubes and result in additional breakage. If a tube breakage occurs, carefully remove the tube holder. Properly dispose of the sample and tube fragments and thoroughly clean both the inside and outside of the tube holder. Insert a new tube cushion (if necessary) and replace the tube holder in the rotor.
4. **Motor and Electrical Maintenance:** The Model 614B uses a brushless A/C motor. It should not need servicing for the life of the centrifuge. The electrical components are selected for high reliability and should not need service.
5. **Keep the Rotor Chamber Clean:** Every six months, or whenever there is a tube breakage, (refer to the note in #3), it may be necessary to remove the rotor and clean the rotor chamber. Follow the instructions on page 5 to remove and re-install the rotor.

CAUTION: Before cleaning, always unplug the line cord from the electrical outlet to eliminate the risk of electric shock.

The rotor chamber, rotor and accessories should be thoroughly cleaned using either isopropyl alcohol, soap and water, or bleach. The use of Fully/Partially Halogenated Hydrocarbons, Ketones, Esters and all other chemicals not prescribed by the manufacturer may cause damage to the rotor and tube holders and should not be used. Apply cleaning solutions with a towel or cloth.

Do not submerge the centrifuge in water or other cleaning solutions as this will cause damage and void your warranty!